



PRECISION SOLUTIONS BY SMC

Enhancing bioprocess efficiency

Expertise
Passion
Automation



SMC solutions for bioreactors

In the evolving landscape of bioprocessing, efficiency, precision, and scalability are crucial for bioreactor manufacturers. At SMC, we specialise in providing high-performance fluid control and automation solutions, ensuring reliable gas, liquid, and temperature control in bioreactor systems. Whether for **seed train development or large-scale production**, our solutions enhance process stability and operational flexibility.

With our expertise in **batch and perfusion systems** as well as **stainless steel (SUS) and glass bioreactor tanks**, we help biopharmaceutical manufacturers meet regulatory demands while optimising performance.

Pinch valves

- For single use line

Thermo control unit (TCU)


- Water temperature control

Digital flow switch

- Monitor the flow of air, Ar, CO₂, N₂, O₂

Flow controller

- Flow control for air, Ar, CO₂, N₂, O₂

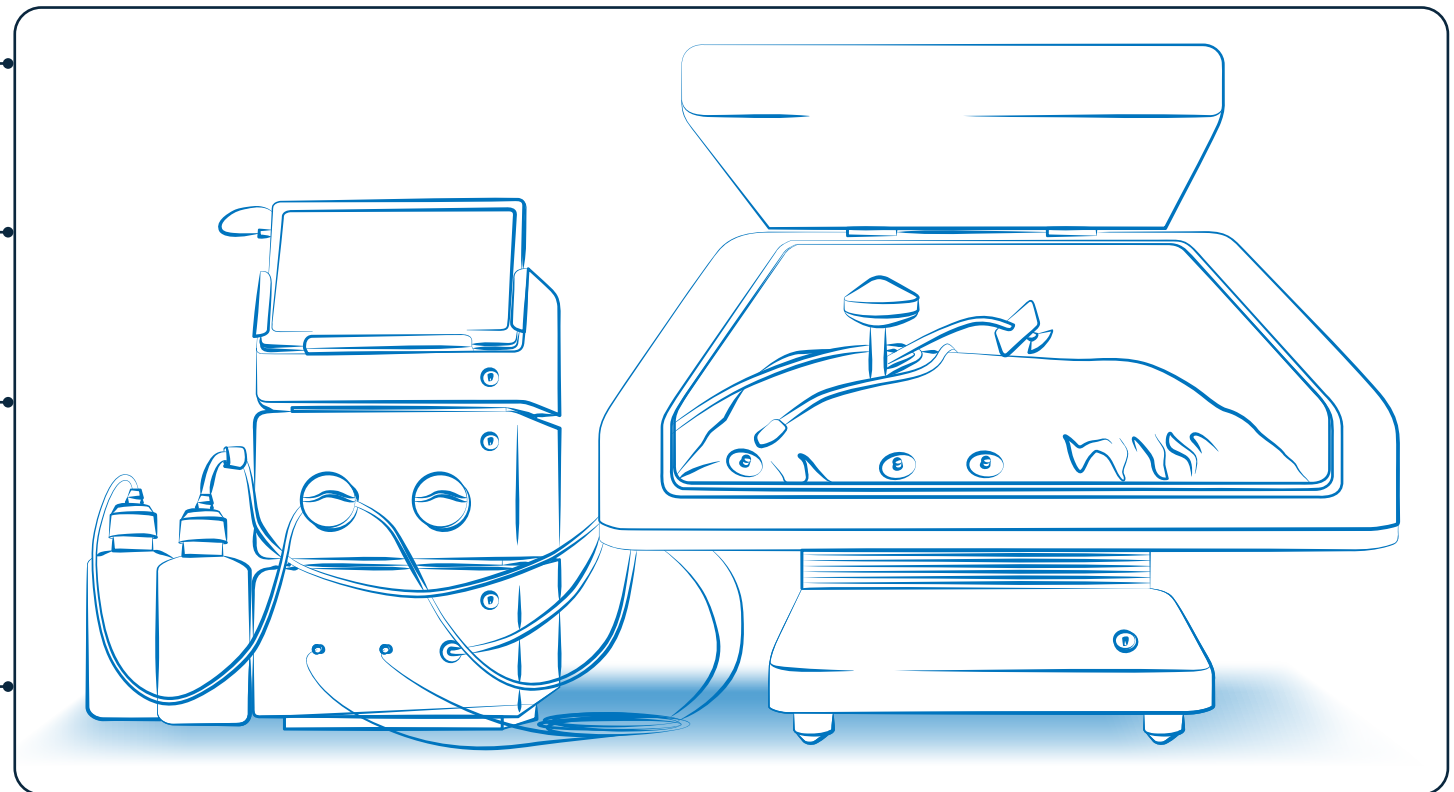

Clean regulator

- Pressure regulation for air, Ar, CO₂, N₂, O₂

Sintered metal element

- For sparger bubbling applications

Directional control valve

- Space saving and improved operability

Bioreactors

Upstream/Seed train

The seed train process, the initial phase of cell culture, is critical for establishing a robust and scalable bioprocess. Bioreactor manufacturers recognise the importance of precise control and reliable components within seed train systems to ensure optimal cell growth and viability. Leveraging its extensive experience in the bioprocessing industry, SMC offers a comprehensive portfolio of solutions designed to meet these demanding requirements



Bio process – Cell culture


Batch/Perfusion bag

SMC offers advanced control solutions for **single-use bioreactors (SUBs)**, ensuring seamless **batch and perfusion processing**. Our solutions improve process reliability while maintaining sterility, making them ideal for **cell culture, fermentation, and biologics production**.

Key benefits:

- **Optimised gas & liquid control** – Precision regulation of **oxygen (O₂), nitrogen (N₂), carbon dioxide (CO₂), and air** for enhanced cell growth
- **Hydrogen technology** – Single-use line compatibility with easy tube replacement for contamination-free operation
- **Temperature stability** – Advanced **thermo-control units (TCUs)** provide precise **cooling and heating** for process consistency
- **Scalability** – Suitable for **R&D, pilot-scale, and large-scale biomanufacturing**.

Our innovative automation solutions ensure that your batch and perfusion processes meet the highest industry standards for both performance and efficiency.

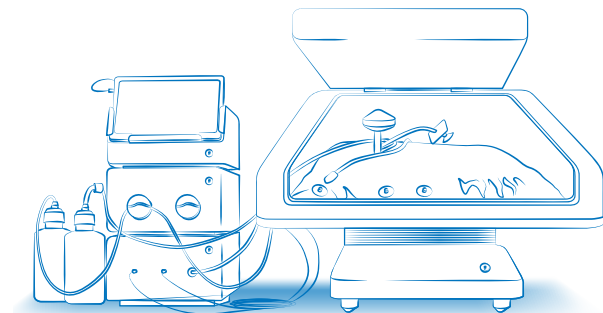


Small & medium size tanks

For bioreactor manufacturers requiring **durable, multi-use or single-use systems**, SMC offers a range of products for stainless steel (SUS), glass and single-use tanks designed for **cell culture, fermentation, and microbial applications**. Our solutions provide precise control over critical bioprocessing parameters, ensuring consistent and repeatable results.

Key benefits:

- **High-purity gas control** – Clean regulators and flow controllers ensure contamination-free sparging of process gases (**O₂, CO₂, N₂, air**)
- **Advanced temperature management** – Thermo-control systems with ± 0.1 °C stability for precise culture conditions
- **Reliable valve solutions** – Low-particle generation **angle seat valves** for seamless integration into stainless steel bioreactor lines
- **Flexible scalability** – Suitable for **pilot-scale to full production-scale applications**, ensuring smooth scale-up.

With SMC's fluid control and automation expertise, bioreactor manufacturers can achieve consistent product quality, process efficiency, and regulatory compliance.

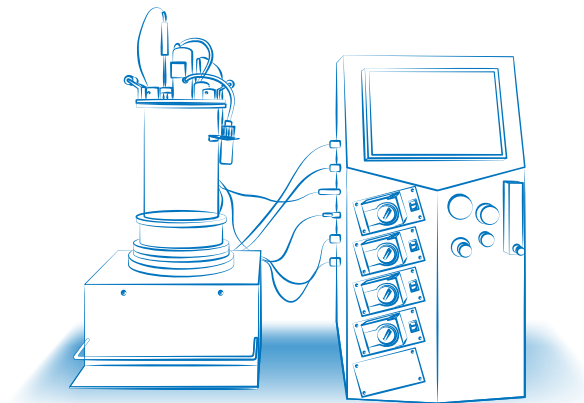

Temperature control

Precise temperature control is essential for maintaining optimal conditions in bioreactors, ensuring cell viability and process efficiency. At SMC, we offer advanced solutions such as high-stability thermal control units, Peltier-based systems, and high-capacity cooling options. Our products provide reliable and precise temperature regulation, meeting the stringent standards of the biopharmaceutical industry. For applications requiring environmental compliance, we offer selected models that are compliant with F-GAS regulations.



HRS Series +

Thermo-control unit

- Up to 38 kw cooling capacity
- Worldwide service location
- Specials available as required.



HRS Series +

Low temp thermo control unit

- Set temperature range and stability
HRS-X158/X164: -10 to 40 °C with ± 0.1 °C.



HEF Series +

Peltier type thermo-control unit

- Low-noise design (37 dB)
- Set temperature range: 10 to 60 °C
- Temperature stability: ± 0.1 °C.



HRSC Series +

Non-F-gas (CO₂ refrigerant) refrigerated thermo-chiller

- Up to 11.5 Kw cooling capacity
- Uses natural refrigerant (CO₂) with GWP=1
- Temperature stability: ± 0.1 °C.

Liquid control

Efficient liquid handling in a bioreactor is key for media feed, collection, and mixing processes. SMC provides a wide range of pinch valves for single-use systems and flow control valves designed for the biotechnology industry. Our solutions ensure precise, contamination-free flow control, adaptable to various applications.



JSX/JSXD/JSXZ/JSXU Series +

2-port valve

- Space-saving and light-weight
- Power consumption: 14 % reduction
- Stainless body and coil.



LVM Series +

Isolated structure valve

- Low particle generation, oil-free, metal-free
- Available with EPDM, FKM, Kalrez®
- Fluid area is separated from solenoid.



JSB Series +

Angle seat valve

- High flow with angle seat design
- Air, water, steam
- Suitable for stainless line with multi-use bioreactor.



LPV Series +

Pinch valve solenoid type


- 2-port and 3-port type
- Easy tube replacement for single-use tube.



KQG2 Series + KFG2 Series +


Stainless fitting

- SUS316 material
- KQG2 one touch type: -5 to 150 °C
- KFG2 Insert type: -65 to 260 °C.

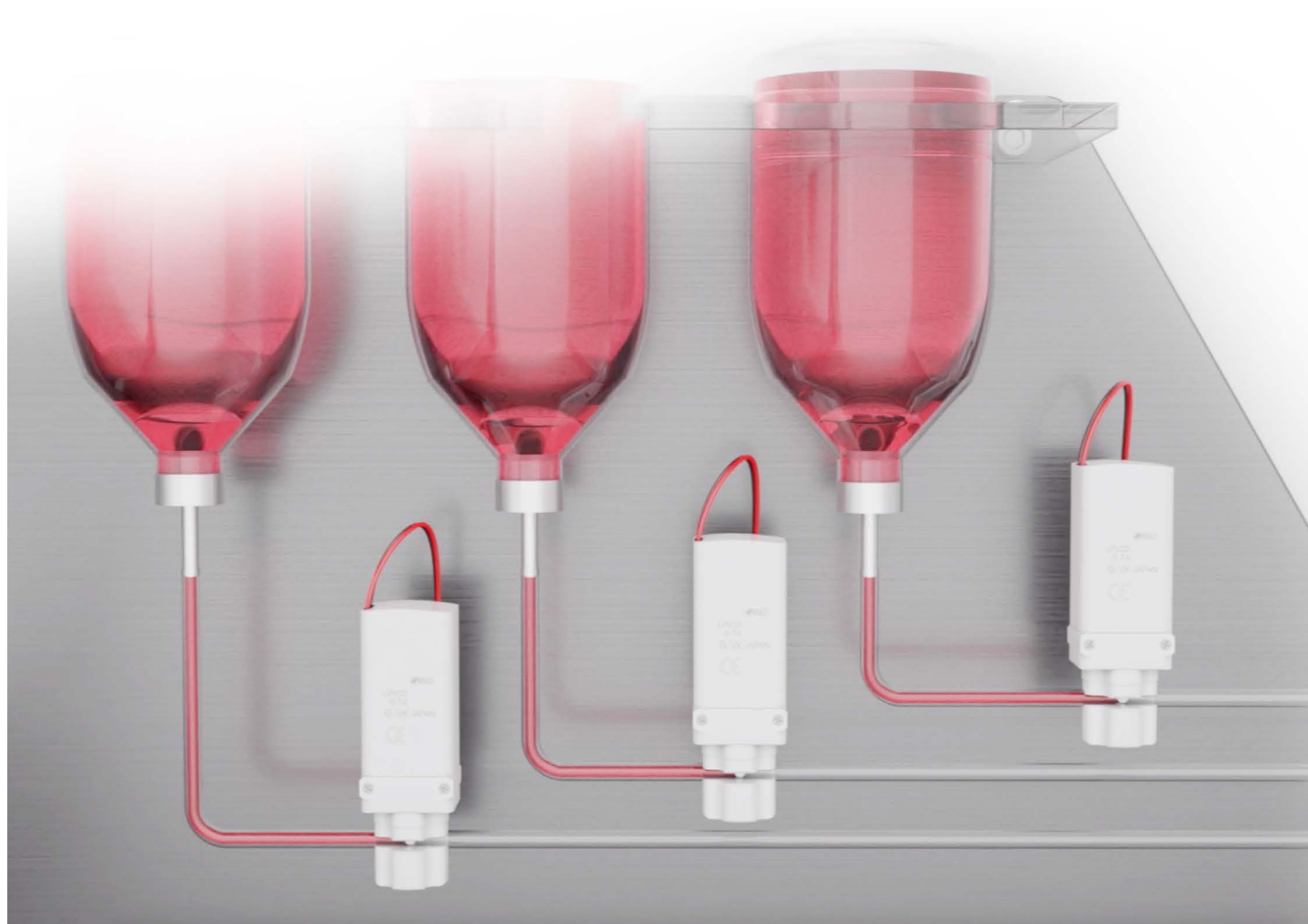

Gas control

Gas supply and control in a bioreactor directly impact oxygenation and cellular metabolic reactions. With our advanced valve technology, pressure regulators, and sintered metal filters, SMC enables precise control of gases such as air, CO₂, O₂, and N₂, optimising oxygen transfer and the efficiency of fermentation and cell culture processes.

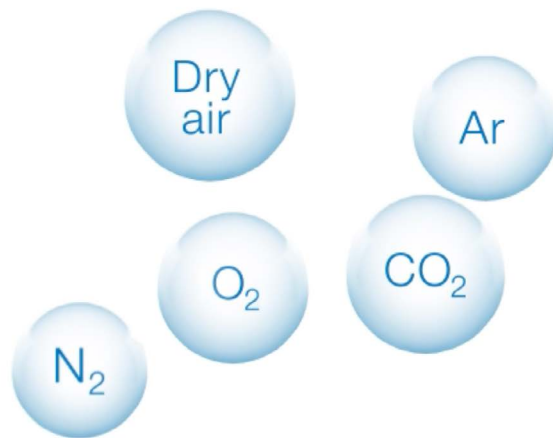

PFC710 Series +

Flow controller

- Flow control for air, Ar, CO₂, N₂, O₂
- Time to flow set: 0.5 seconds or less
- Repeatability: ± 1 % F.S.



PF2M Series +

Digital flow switch

- 2 to 200 l/min range
- Air, Ar, CO₂, N₂, O₂
- Grease-free.



SR□ Series +

Clean regulator

- Pressure regulation for air, Ar, CO₂, N₂, O₂
- Contamination-controlled SUS steel
- Oil-free.



ES□ Series +

Sintered metal filter element

- Space-saving and light weight
- Stainless material
- For sparger bubbling applications.

Device communication

The seamless integration and automation of bioreactor systems require reliable communication solutions. SMC offers serial interfaces with protocols such as Ethernet/IP™, PROFINET, and IO-Link, ensuring real-time monitoring and connectivity. Our solutions enable smooth integration with automation systems, enhancing process control and traceability.



EXW1/EX600-W Series +

Serial interface

- Wireless system option
- Various protocols available (Ethernet/IP™, PROFINET, IO-Link, etc).

IO-Link

EtherNet/IP™

EtherCAT®

PROFINET®
BUS

PROFINET®

DeviceNet™

CC-Link V2



Air line components

The quality of compressed air is a critical factor in bioreactor systems, directly affecting sterility and process performance. SMC provides high-performance pneumatic regulators, air treatment units, and control valves designed to maintain a clean and stable air supply. Our solutions ensure operational reliability and energy efficiency in critical applications.



AC-D Series +

Air combination units for air supply

- Urban white, a clean look for your equipment
- Modular design
- Element and bowl form one piece.



SY Series +

Directional control valve


- Space-saving and improved operability
- Power consumption of 0.1 W.



JSY Series +

Directional individual valve

- Flow size reduction
- Dual 3-port valves available for supplying pilot air
- Low power consumption (0.1 W).

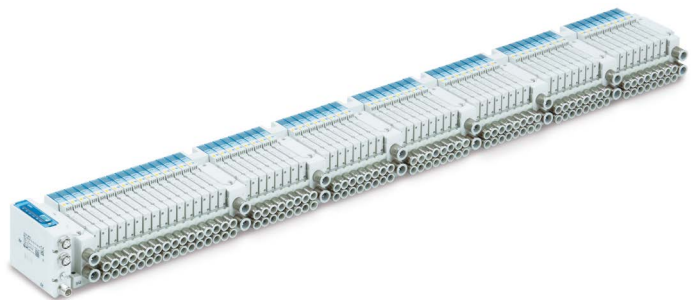

JSY3000-L Series +

64-station compatible manifold


- Valve stations: from 4 to 64
- Reduced number of SI units
- Reduced wiring and wiring work.



ITV Series +

Electro-pneumatic regulator

- Stepless air pressure control
- Linearity: ± 1 % or less (F.S.)
- Hysteresis: 0.5 % or less (F.S.).

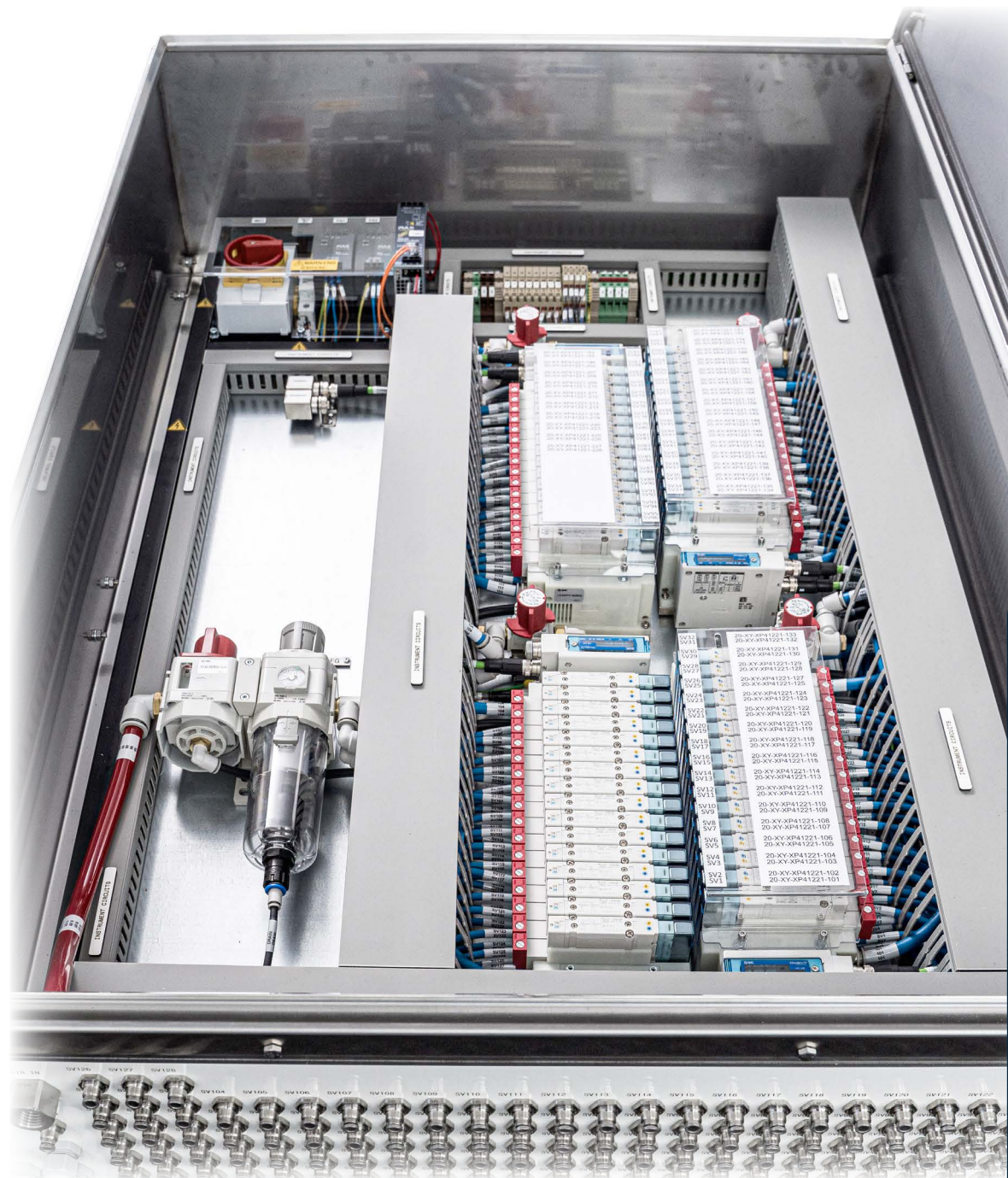


Our support network

SMC's worldwide commitment

One of the things we do best at SMC is **being close to our customers**. Local support, on a global scale.

With **support** in over **500 locations** across **80 countries** and regions **worldwide**, our sales force of **7000 experts** maintains **close communication with customers**.

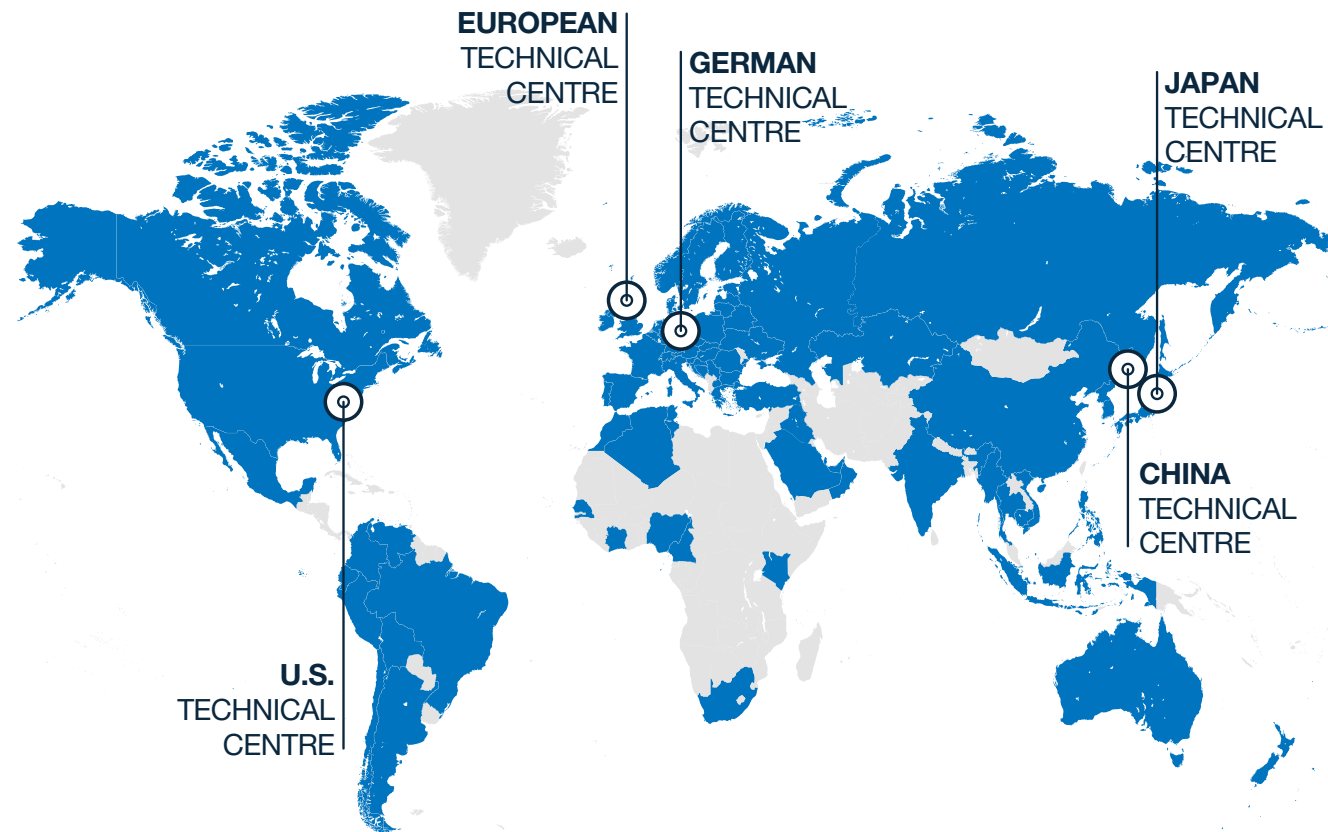

SMC Business Continuity Plan

Sustainable growth also means ensuring uninterrupted operations

We are committed to ensuring that SMC is prepared for any emergency and that our business activities will not stop in the event of such circumstances. SMC aims to fulfil our product supply responsibilities and maintain our customers' trust by contributing to both sustainable growth and the expansion of technological innovations.

SMC, as a comprehensive manufacturer of automatic control equipment that supports automation, is able to promptly provide products that meet our customers' needs anywhere in the world.

Production BCP

Ensure customer order fulfilment

Reliable delivery for you thanks to our 9 global logistic centres and 38 production sites worldwide, 10 of which are located in Europe. Moreover, flexibility to rapidly respond to any sudden change in the manufacturing environment.

Finance BCP

Safe & Solid financial base

In the event of an emergency, SMC can provide a safe and solid financial base (with cash, deposits, and equity capital) that will sufficiently cover the working capital and funds needed to rebuild buildings and the equipment required for business continuity. This is done to provide peace of mind to our customers and workers alike.

Information security BCP

Vital data kept safe

Strengthen information security for protection against computer viruses and cyberattacks, plus the installation of data centres to establish a disaster recovery system. Your information is safe with us.

Engineering BCP

Consistent technical support

2,100 engineers at our 5 technical centres around the globe (2 in Europe – Germany and UK).

Sales BCP

Consistent sales support

7,000 sales engineers worldwide ready to recommend the best solution for you. Over 80 global locations to make sure that wherever you are, we are there too.



Discover more



SMC Corporation

1-5-5, Kyobashi,
Chuo-ku, Tokyo
104-0031, Japan
Telephone: 03-6628-3000
<https://www.smcworld.com>

Austria	+43 (0)2262622800	www.smc.at	office.at@smc.com
Belgium	+32 (0)33551464	www.smc.be	info@smc.be
Bulgaria	+359 (0)2807670	www.smc.bg	sales.bg@smc.com
Croatia	+385 (0)13707288	www.smc.hr	sales.hr@smc.com
Czech Republic	+420 541424611	www.smc.cz	office.at@smc.com
Denmark	+45 70252900	www.smc.dk.com	smc.dk@smc.com
Estonia	+372 651 0370	www.smc.ee	info.ee@smc.com
Finland	+358 207513513	www.smc.fi	smc.fi@smc.com
France	+33 (0)164761000	www.smc-france.fr	supportclient.fr@smc.com
Germany	+49 (0)61034020	www.smc.de	info.de@smc.com

Greece	+30 210 2717265	www.smcHELLAS.gr	sales@smcHELLAS.gr
Hungary	+36 23513000	www.smc.hu	office.hu@smc.com
Ireland	+353 (0)14039000	www.smcautomation.ie	technical.ie@smc.com
Italy	+39 03990691	www.smcitalia.it	mailbox.it@smc.com
Latvia	+371 67817700	www.smc.lv	info.lv@smc.com
Lithuania	+370 5 2308118	www.smc.lt	info.lt@smc.com
Netherlands	+31 (0)205318888	www.smc.nl	info@smc.nl
Norway	+47 67129020	www.smc-norge.no	post.no@smc.com
Poland	+48 22 344 40 00	www.smc.pl	office.pl@smc.com
Portugal	+351 214724500	www.smc.eu	apoiocliente.pt@smc.com

Romania	+40 213205111	www.smcromania.ro	office.ro@smc.com
Russia	+7 (812)3036600	www.smc.eu	sales@smcru.com
Slovakia	+421 (0)413213212	www.smc.sk	sales.sk@smc.com
Slovenia	+386 (0)73885412	www.smc.si	office.si@smc.com
Spain	+34 945184100	www.smc.eu	post.es@smc.com
Sweden	+46 (0)86031240	www.smc.nu	order.se@smc.com
Switzerland	+41 (0)523963131	www.smc.ch	helpcenter.ch@smc.com
Turkey	+90 212 489 0 440	www.smc-turkey.com.tr	satis.tr@smc.com
UK	+44 (0)845 121 5122	www.smc.uk	sales.gb@smc.com
South Africa	+27 10 900 1233	www.smcza.co.za	Sales.za@smc.com

www.smc.eu

Release DW
BIOREACTOR-A-EN

SPECIFICATIONS ARE SUBJECT TO CHANGE WITHOUT PRIOR NOTICE AND ANY OBLIGATION ON THE PART OF THE MANUFACTURER